


Vitreotomy Cutters - 23Ga

Posterior Vitrectomy Cutters, 23Ga

5172 30 PSI	Alcon: Accurus - 2500 cpm DORC: Associate - 2500 cpm
5216 25 PSI	DORC: Harmony Total TTC - 1500 cpm
5171 25 PSI	B+L: Millennium, Premiere - 750 cpm Carl Zeiss Meditec/IOL Tech: Pentasys - 1500cpm Carl Zeiss Meditec: Visalis 500 Fritz Ruck: Pentasys, Pentasys 2 - 2500 cpm Croma Pharma/Corneal: Open Phaco - 1500 cpm Optikon: Pulsar ² , Assistant - 2000 cpm Nidek: CV-12000, CV-6000, CV-7000 - 800 cpm Nidek: CV-24000 - 2500 cpm

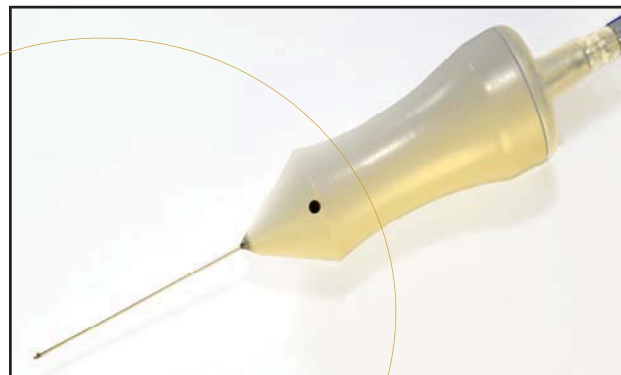

Anterior Vitrectomy Cutters, 23Ga

5208 20 PSI	Alcon: Legacy, Universal AMO: Sovereign, Signature, AMO Plus, Prestige Staar Surgical: Wave Nidek: CV-12000, CV-6000, CV-7000
-----------------------	--

Fiber Optic Probes, 23Ga

5177	Fiber Optic Probe, 23Ga - Alcon Accurus - 6 per box
5178	Wide Angle Fiber Optic Probe, 23Ga - Alcon Accurus - 6 per box
5253	Fiber Optic Probe, 23Ga - Alcon Constellation - 6 per box
5254	Wide Angle Fiber Optic Probe, 23Ga - Alcon Constellation - 6 per box
5179	Fiber Optic Probe, 23Ga - B+L Millennium - 6 per box
5180	Wide Angle Fiber Optic Probe, 23Ga - B+L Millennium - 6 per box
5247	Fiber Optic Probe, 23Ga - B+L Stellaris PC - 6 per box
5237	Wide Angle Fiber Optic Probe, 23Ga - B+L Stellaris PC - 6 per box
5163	Fiber Optic Probe, 23Ga - Ocutome connector - 6 per box
5169	Wide Angle Fiber Optic Probe, 23Ga - Ocutome connector - 6 per box
5210	Shielded Fiber Optic Probe, 23Ga - Ocutome connector - 6 per box
5249	Chandelier Fiber with Adjustable Ring, 25Ga/0.5mm - B+L Stellaris PC - 6 per box
5250	Chandelier Fiber with Adjustable Ring, 25Ga/0.5mm - Alcon Constellation - 6 per box
5228	Chandelier Fiber with Adjustable Ring, 25Ga/0.5mm - Ocutome 5mm connector - 6 per box

One-Step Trocar System - 23Ga

Available January 2012

- Self-sealing wound construction
- Ultra sharp knife for smooth incision
- Unique closure valve system
- 3.0 + 4.0 mm scleral marker

5229 One-Step Trocar System, 23Ga


Each package includes:

- Three Ultra sharp Trocar Cannula Knives
- Three Preloaded self-sealing Trocar Cannulas
- One Universal Infusion Line




Trocar Cannula Knife

Knife Inserter consist of ultra sharp blade with preloaded Self-Sealing Trocar Cannula

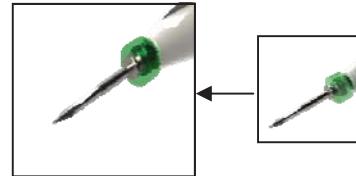

Self-Sealing Trocar Cannula

Stainless steel Trocar Cannula is complete with self-sealing silicone valve, eliminating need for plugs. Valve offers superior sealing and smooth instruments extraction

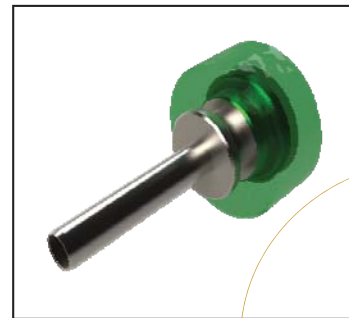

Universal Infusion Line

- Instantly connects to Trocar Cannula
- Designed to offer maximum infusion flow



Two-Step Trocar System - 23Ga

Available January 2012

Second Generation

- Self-sealing wound construction
- Ultra sharp knife for smooth incision
- Unique closure valve system
- 3.0 + 4.0 mm. scleral marker

5243 Two-Step Trocar System, 23Ga

Each package includes:

- Three Trocar Cannula Inserters
- Three Preloaded self-sealing Trocar Cannulas
- One Trocar Knife
- One Universal Infusion Line




Trocar Cannula Inserter

Inserter assembled with preloaded Self-sealing Trocar Cannula



Self-Sealing Trocar Cannula

Stainless steel Trocar Cannula is complete with self-sealing silicone valve, eliminating need for plugs. Valve offers superior sealing and smooth instruments extraction

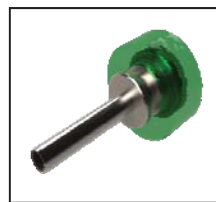

Trocar Knife

Trocar Knife consists of ultra sharp stiletto mounted to plastic handle



Universal Infusion Line

- Instantly connects to Trocar Cannula
- Designed to offer maximum infusion flow



Two-Step Trocar System - 23Ga

First Generation

Trocar Cannula Inserter works together with Trocar Knife to accomplish Self-Sealing Trocar Cannula placement without using fixation plate.



5181 Two Step Trocar System, 23Ga


Each package includes:

- Three Trocar Cannula Inserters
- Preloaded Self-sealing Trocar Cannulas
- One Trocar Knife
- One Infusion Line

Watch video clip on [MicroVision website](#)

Trocar Cannula Inserter

Trocar Cannula Inserter comes with Cannula pre-installed. Delivered in retracted position.

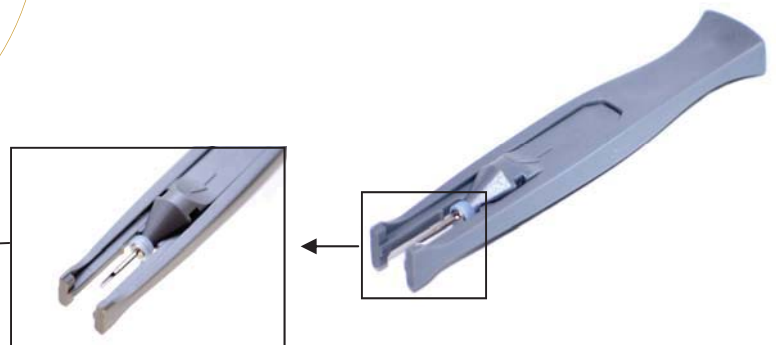


Self-Sealing Trocar Cannula

Stainless steel Trocar Cannula is complete with self-sealing cap, eliminating need for plugs.



Trocar Knife

Trocar Knife consists of sharp trocar mounted to plastic handle. Delivered in retracted position.

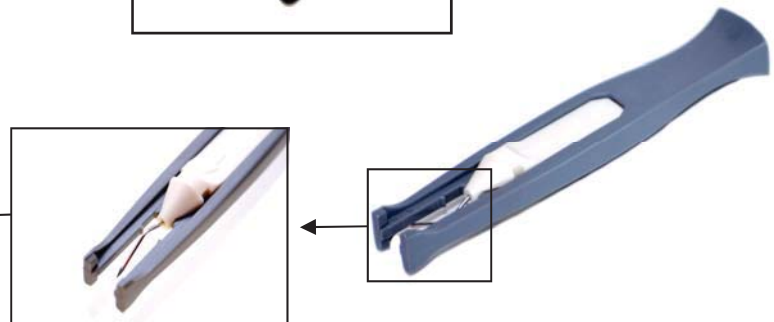

Two-Step Trocar System - 23Ga

Reusable Two Step Trocar System with closure valves, 23Ga

- 5173-1 Package includes:
- Trocar Cannula with closure valves - 3 pcs
 - Universal Infusion Line - 1 pc
 - Blunt Cannula Inserter - 3 pcs
 - Fixation Plate - 1 pc
 - Loading Forceps - 1 pc




Reusable Two Step Trocar System with closure plugs, 23Ga

- 5173 Package includes:
- Trocar Cannula - 3 pcs
 - Trocar Cannula Plug - 3 pcs
 - Universal Infusion Line - 1 pc
 - Blunt Cannula Inserter - 3 pcs
 - Fixation Plate - 1 pc
 - Loading Forceps - 1 pc



[Watch video clip on MicroVision website](#)


5186 Loading Forceps

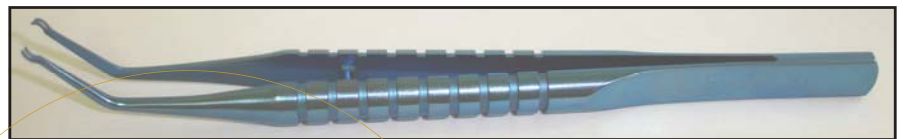

5187 Blunt Cannula Inserter



5188 Fixation Plate

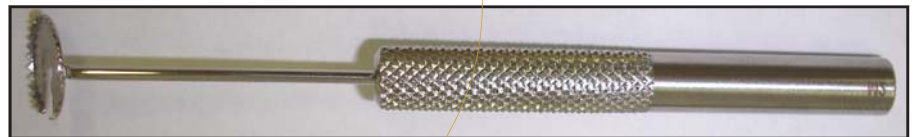

5202 Universal Infusion Line
Packaged sterile, 5 per box



Vitrectomy Accessories - 23Ga

Trocar Cannula Set, 23Ga

5190

Each Set includes:

- Instrument Cannula - 2ea
 - Cannula Plug - 2ea
- Sterile, five sets per box



Closure Valves for 23Ga Trocar Cannula

5209

Replacement Closure Valves - 1 set of 4pcs, sterile



5184 Angled MVR Blade, 23Ga - 6 per box




Self-Retaining Sutureless Cannulas, 23Ga

5144 Self-Retaining Sutureless Infusion Cannula 23Ga, 4.0mm

5221 Self-Retaining Sutureless Infusion Cannula 23Ga, 6.0mm

5165 Self-Retaining Sutureless Silicone Oil Cannula 23Ga, 4.0mm

5222 Self-Retaining Sutureless Silicone Oil Cannula 23Ga, 6.0mm

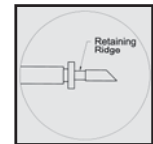

Dual Bore Subretinal BSS Injection Needle, 20Ga/23Ga



- 5257 - Dual Bore BSS Injection Needle, 23Ga with 0.1mm/41Ga tip - 1 per box
 - Enables to control subretinal injection of BSS
 - Dual bore cannula combined with aspiration capability

Viscous Fluid Injection Cannula, 23Ga



- 5248 - Viscous Fluid Injection Cannula 23Ga, 10mm polyimide tip - 5 per box
- Allowing injection of viscous fluids such as silicone oil through 23Ga trocar cannula

Dual Bore PFC Cannula, 23Ga

5203 Dual Bore PFC Cannula, 23Ga - 5 per box

- Simultaneous infusion of heavy liquids and aspiration of intraocular fluids
- Dual bore enables to control constant intraocular pressure during injection

